

Outline


This compact NV-300 requires no inspector and achieves both 5 m far vision and 40 cm near vision measurements.

Features

- 1) Built-in start sensor automatically measures the visual acuity without inspector.
- 2) Voice guide achieves the easy operation by any person.
- 3) Far vision and near vision can be switched automatically. The corrected vision or naked vision is automatically inputted.
- 4) A single chart presented on the center of the visual field can be identified easily and achieves the correct measurement.
- 5) An indication lamp makes it possible to know the general visual acuity from outside.
- 6) High-speed printer is built in.

Specifications

- 1) Visual acuity chart : 4 types, top, bottom, right and left, of Landolt's ring
- 2) Visual acuity value : 12 step; 0.1, 0.2, 0.3, 0.4, 0.5, 0.6, 0.7, 0.8, 0.9, 1.0, 1.2 and 1.5
- 3) Chart distance : 5 m
- 4) Inspection item : 5 m far vision, 40 cm near vision, right, left and both eyes
- 5) Power supply : AC100 / 120 / 230 V 50 / 60 Hz 60 VA
- 6) Lamp : White LED (replacement-free)
- 7) Data output : RS-232C, D-SUB 9 pin
- 8) Dimensions & weight : 220 (W) x 410 (D) x 340 (H) mm, 10 kg

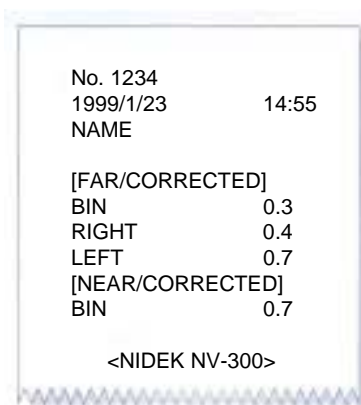

Specifications and design are subject to change without notice for improvement.



NIDEK CO., LTD.

HEAD OFFICE	:34-14, Maehama, Hiroishi Gamagori, Aichi 443-0038, Japan Telephone :81-533-67-6611 Facsimile :81-533-67-6610 URL :http://www.nidek.co.jp
TOKYO OFFICE (International Div.)	:6F Takahashi Bldg., 3-2 Kanda-Jinboucho Chiyoda, Tokyo 101-0051, Japan Telephone :81-3-3288-0571 Facsimile :81-3-3288-0570 URL :http://www.nidek.com
NIDEK INC.	:47651 Washington Drive Fremont, CA 94539, U.S.A. Telephone :1-510-226-5700 :1-800-223-9044 (US only) Facsimile :1-510-226-5750 URL :http://www.usa.nidek.com
NIDEK TECHNOLOGIES AMERICA INC.	:5500 West Friendly Ave. Suite 101 Greensboro, NC 27410, U.S.A. Telephone :1-336-851-0225 :1-888-382-5064 (US only) Facsimile :1-336-851-0917 URL :http://www.nidektech.com
NIDEK SOCIÉTÉ ANONYME	:Europarc 13, rue Auguste Perret 94042 Créteil, France Telephone :33-1-49 80 97 97 Facsimile :33-1-49 80 32 08 URL :http://www.nidek.fr
NIDEK TECHNOLOGIES SRL.	:Via dell'Artigianato, 6/A 35020 Albignasego (Padova), Italy Telephone :39 049 8629200 / 8626399 Facsimile :39 049 8626824 URL :http://www.nidektechnologies.it